


## KEN IWD 521

### Washer-disinfector for CSSD-departments

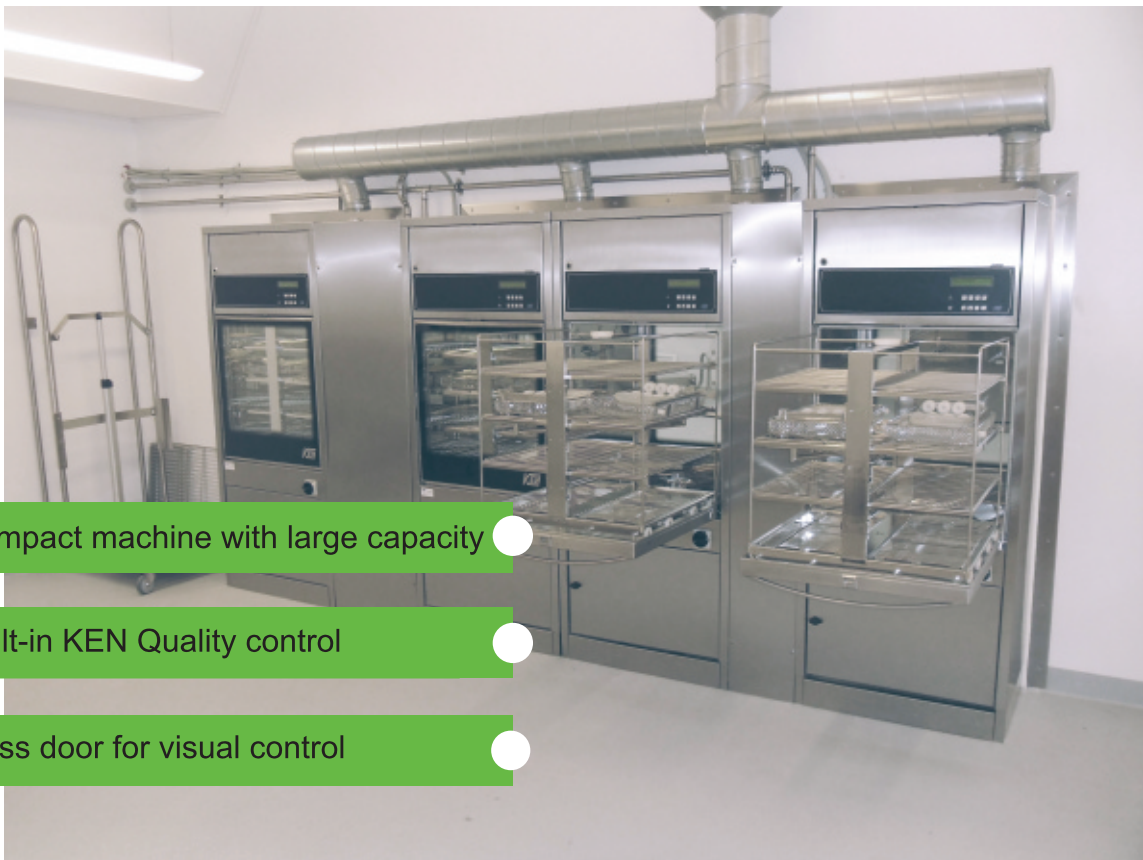

Compact machine with large capacity

Built-in KEN Quality control

Glass door for visual control

For more information consult - [www.ken.dk](http://www.ken.dk)

### The flexible machine with considerable capacity



KEN IWD 521 is particularly suitable for use in sterile centres as well as in surgical and anaesthetic wards where there is a high quality, capacity, flexibility and full process documentation requirement.

KEN IWD 521 can wash at up to 5 levels and is available both as a pass-through machine and with one door.

The pass-through model ensures effective separation of an unclean and clean side thus optimizing the level of hygiene and reducing the risk of hospital infections.

The machine complies with EN ISO 15883.



## Service

The machine has built-in KEN Quality control for optimum reliability. It demands a quality assurance check every 24 hours. The construction is service-friendly with easy access to all components.

All KEN machines conform to present EU directives and all hospital machines have been constructed on the basis of EN ISO 15883 including for example the testing of its cleaning capacity according to SIS-TR 3:2001 (formerly SPRI).



## Advantages

- Computerized monitoring for always correct inside and outside washing and disinfection
- Optional KEN NetCom or XML Connect software solutions
- Washing of instruments at up to 6 levels at the same time
- 12 DIN baskets
- User-friendly and ergonomically correct front operation
- Front display in plain language
- Low noise level
- Low consumption of energy and chemicals - environmentally friendly and economical
- Minimum of maintenance and servicing
- Possibility of re-programming on the spot as well as the addition of new programmes
- Possible connection to built-in printer for optimum process documentation

## Accessories

Above are examples of accessories/equipment for this machine. Look for equipment with the aid of the equipment selector at [www.ken.dk](http://www.ken.dk) or contact KEN for further information including customized special solutions

## Contact a KEN specialist

KEN will gladly help you find the right solution. Consult [www.ken.dk](http://www.ken.dk) to find your KEN dealer.

Tel: +45 62 63 10 91  
E-mail: [ken@ken.dk](mailto:ken@ken.dk)  
Web: [www.ken.dk](http://www.ken.dk)

## Technical specifications

|                                      |   |
|--------------------------------------|---|
| Number of washing levels             | 5-6   |
| Capacity, baskets                    | 10-12 DIN-baskets 480 x 250 mm.                       |
| Volume L                             | 270 L   |
| Programme time, standard programme   | 30-55 min. excl. drying                               |
| Power connection                     | 3 x 400 V, N, E, 50Hz                                 |
| Connected load / initial fuse        | 13,6kW / 25A  |
| Water connection                     | Cold untreated water: 1 pcs. 3/4" RG, 5-30°C, 0-30°dH |
| Water connection                     | Hot soft water: 1 pcs. 3/4" RG, 0-3°dH                |
| Water connection                     | Demin. water: (Extra) 1 pcs. 3/4" RG 5-70°C           |
| Outlet                               | Ø 50mm.   |
| Water pressure                       | 200-800kPa  |
| Amount of water (approx.) per intake | 20 L  |
| Measurements machine W x D x H       | 746mm (condensator+50mm)x740x1820 mm                  |
| Dosage of instrument lubricant/acid  | Extra   |
| Chamber size W x D x H               | 640 x 640 x 670 mm                                    |
| Pump capacity                        | 600 L/min.  |

**We make a difference - locally and globally**